

oerlikon
leybold vacuum

UNIVEX

High Vacuum Experimentation Systeme

Vacuum Coating Technology
Custom Systems

184.01.02

Excerpt from the Oerlikon Leybold Vacuum Full Line Catalog

Product Section C19

Edition May 2007

Contents

General

General	C19.02
Applications and Accessories / Process Components	C19.03

Products

Bell Jar System	
UNIVEX 300	C19.04
Options for the UNIVEX 300	C19.07
Door System	
UNIVEX 350	C19.08
UNIVEX 350 G	C19.12
UNIVEX 450 B	C19.16
Special Systems	
UNIVEX 450 for Dactyloscopy	C19.20
UNIVEX 450 C / Test Systems with a Vacuum Chamber	C19.21

Accessories

Standard Accessories for UNIVEX Systems	C19.22
Components	
for Thermal Evaporation of High Melting Point Materials (metals)	C19.24
for Thermal Evaporation of Low Melting Point (organic) Materials	C19.26
for Electron-Beam Evaporation	C19.26
for Sputtering	C19.27
for Glow Discharge Cleaning	C19.28
for Film Thickness Measurements	C19.29
Substrate Manipulation Solutions	C19.30

General

General

The UNIVEX multi-purpose experimentation systems were developed by Oerlikon Leybold Vacuum for applications in research and development, as well as for setting up pilot production systems.

Their range of applications focuses chiefly on vacuum coating technology as well as vacuum process engineering experiments.

The universal experimentation systems from Oerlikon Leybold Vacuum are

based on a modern modular concept suited for customised expansion. The high vacuum pumps are fitted horizontally at the side of the base panel, respectively the vacuum chamber.

Application and Accessories

	UNIVEX 300	UNIVEX 350	UNIVEX 350 G	UNIVEX 450 B	Custom versions
Applications	Bell jar system		Door system		
Passive components	■	■	■	■	■
Sensor technology	■	■	■	■	■
Opto-electronics	■	■	■	■	■
Tribology	■	■	■	■	■
Soldering	■	■	■	■	■
Dactyloscopy					■
Glove box applications	■	■	■	■	■
Special applications	■	■	■	■	■
Thermal conduction experiments	■				■

Accessories / Process Components

	Page	UNIVEX 300	UNIVEX 350	UNIVEX 350 G	UNIVEX 450 B	Custom versions
Standard accessories	C19.22	■	■	■	■	■
Thermal evaporation	C19.24	■	■	■	■	■
Electron-beam evaporation	C19.26	■	■	■	■	■
DC sputtering	C19.27	■	■	■	■	■
RF sputtering	C19.27		■	■	■	■
Glow discharge cleaning	C19.28	■	■	■	■	■
Film thickness measurement	C19.29	■	■	■	■	■
Substrate manipulation (rotation, heating, cooling, bias)	C19.30	■	■	■	■	■

Products

UNIVEX 300



UNIVEX 300, typical arrangement with stainless steel bell jar and process components

Bench System with a Vacuum Bell Jar (300 mm diameter)

Typical Applications

- Vacuum coating in research and development
- Special experiments

Basic Unit

- The pump system and the electrical supply system are housed in a 19" rack cabinet
- Moreover, the 19" cabinet provides space for a vacuum gauge and a thickness measuring instrument as well as power supply units for the process components
- The basic unit may be placed on a bench top

Vacuum Chamber

- The base plate is attached to the lateral intake port of the basic unit
- Either a vacuum chamber made of stainless steel or glass may be placed on the base plate

Pump System

- The standard pumping equipment comprises a TRIVAC D 8 B two-stage rotary vane pump and a TURBOVAC 361 turbomolecular pump

Vacuum Measurement

- Depending on the type of application, a combination vacuum gauge operating according to the cold cathode or hot cathode principle may be installed

Advantages to the User


- Modular system design
- Any kind of process component may be installed

- Process components may be retrofitted without problems
- Free access to vacuum bell jar, base plate and chamber units


Options

- Upon request, the UNIVEX 300 may be equipped with a manual lifting facility for the bell jar or also a PLC controller for process automation
- For processes in which aggressive media need to be pumped, a seal gas version of the turbomolecular pump and a rotary vane pump with a special oil filling can be supplied
- For particularly sensitive processes, a dry compressing backing pump like the SCROLLVAC SC 30 D can be used

More information on these options is available upon request.

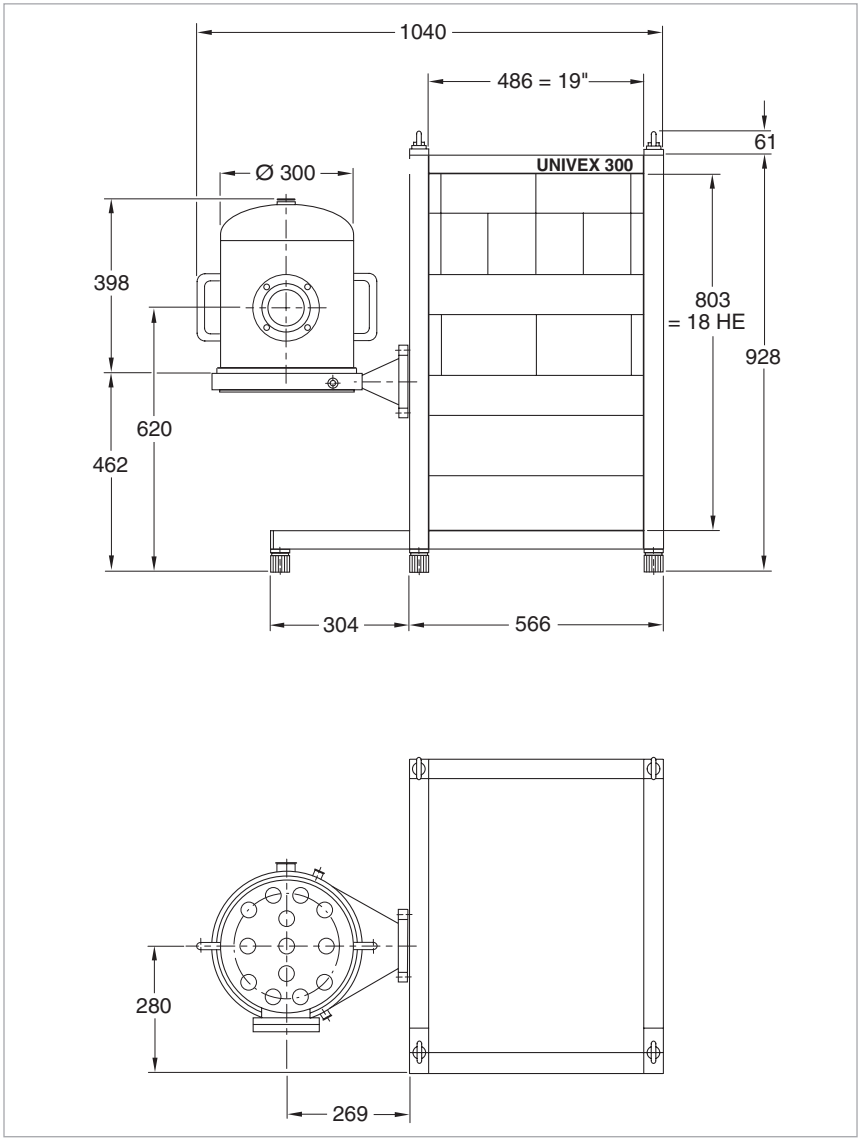

Dimensional drawing for the UNIVEX 300 with base plate and stainless steel bell jar

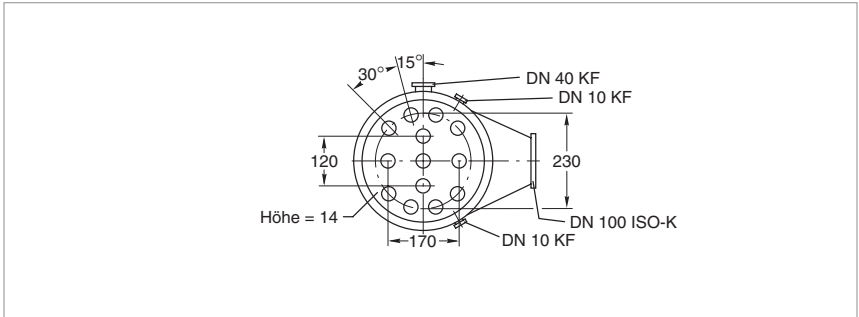

Dimensional drawing for the base plate on the UNIVEX 300

Technical Data

UNIVEX 300

Stainless Steel with Viewing Window (optional: Pyrex glass)

Vacuum chamber		Base panel with bell jar fitted from the top
Diameter	mm	300
Base plate		
Material		stainless steel
High vacuum connection flange	DN	100 ISO-K, attached at the side
Dimensions (H x dia.)	mm	60 x 350
Installation bores	mm	∅ 34,5 (13 x) ¹⁾
Further side flanges	DN	2 x 10 ISO-KF, 1 x 40 ISO-KF ¹⁾
Vacuum bell jar		Stainless steel with viewing window (optional: Pyrex glass)
High vacuum pump		TURBOVAC 361
Nominal pumping speed for N ₂	l x s ⁻¹	345
Power supply		TURBOTRONIK NT 20
Backing pump		TRIVAC D 8 B
Nominal pumping speed	m ³ x h ⁻¹	9.7
Controller		Power supply with main switch slot module
Electrical connection		230 V, 50 Hz, max. 16 A ^{2), 3)}
Cooling water connection		
Hose DN 10	bar	4 to 7
Cooling water consumption, approx.		
l x min ⁻¹		1 ³⁾
Weight, approx.	kg	170 ³⁾

Ordering Information

UNIVEX 300

UNIVEX 300	upon request
------------	--------------


¹⁾ Standard configuration, other hole patterns/flanges upon request

²⁾ Other voltages and frequencies upon request

³⁾ Without chamber installations/process components

Options for the UNIVEX 300

Pyrex Glas Bell Jar (Vacuum Chamber)

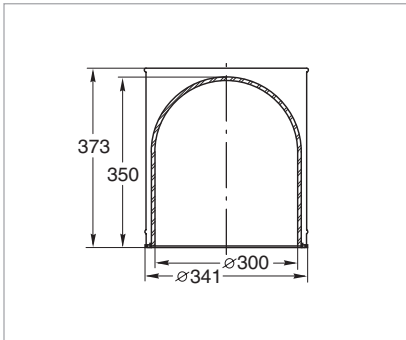

Dimensional drawing for the glass bell jar with implosion protection

Technical Data

Dimensions (H x dia.)	mm	350 x 300
Height, cylindrical section	mm	200
Seal		FPM
Weight	kg	5.6

Bell Jar, Pyrex Glas


Ordering Information

Bell jar, Pyrex glas	upon request ¹⁾
----------------------	-----------------------------------

Bell Jar, Pyrex Glas

¹⁾ With punched steel cover for implosion protection

Stainless Steel Bell Jar (Vacuum Chamber)

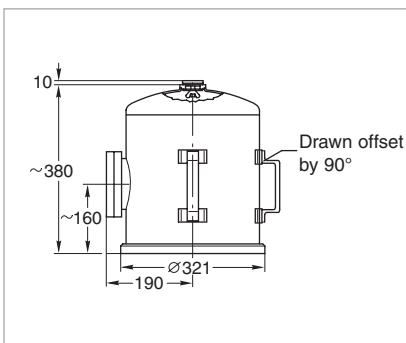

Dimensional drawing for the stainless steel bell jar

Technical Data

Dimensions (H x dia.)	mm	380 x 300
Height, cylindrical section	mm	300
Seal		FPM
Weight	kg	9.6

Bell Jar, Stainless Steel

Ordering Information

Bell jar, stainless steel	upon request ²⁾
---------------------------	-----------------------------------

Bell Jar, Stainless Steel

²⁾ With DN 100 viewing window and 2 carrying handles; hole at the top (34.5 mm dia.)

UNIVEX 350



UNIVEX 350



UNIVEX 350, custom version with higher vacuum chamber and on castors



UNIVEX 350, custom version with second coating module

Laboratory System with Cubic Vacuum Chamber and Front Door (350 mm diameter)

Typical Applications

- Vacuum coating in research and development
- Pre-production trials
- Special experiments

Design

- The UNIVEX 350 consists of two separable 19" rack mount cabinets
- The process chamber and the pump system are accommodated in one cabinet
- Accommodated in the second cabinet section is the electric power supply with the PLC based system controller with graphic touchscreen. Moreover, the supply units for the process components are also accommodated in this cabinet section

Vacuum Chamber

- The base plate is attached to the base frame
- The door is equipped with a viewing window
- Bottom plate and lid are provided with installation holes
- Additional flanges at the side for pump system and process components

Pump System

- The standard pumping equipment comprises a TRIVAC D 25 B two-stage rotary vane pump and a TURBOVAC 600 C turbomolecular pump

Vacuum Measurement


- Depending on the type of application, a combination vacuum gauge operating according to the cold cathode or hot cathode principle may be installed

Advantages to the User


- Modular system design
- Freedom of installation and retrofitting of process components without problems
- Vacuum chamber with front door for highly flexible use
- Convenient access to the chamber installations
- Very simple to operate and use via programmable control
- For installation into clean-room walls
- For RF sputtering
- Pump system adapted to the individual process

Options

- A water-cooled vacuum chamber can also be supplied
- Evaporation protection plates which may be easily disassembled are available
- A water-cooled vacuum chamber can also be supplied
- Evaporation protection plates which may be easily disassembled are available
- Fitting of a second coating module with a vacuum chamber (fitted to the right of the electrical cabinet) is provided for
- For processes which develop increased quantities of gas or which require low operating pressures, the UNIVEX 350 can also be equipped with a turbomolecular pump having a higher pumping speed (TURBOVAC 1000, for example) or with cryopumps
- For processes which involve pumping of aggressive media, a barrier gas version of the turbomolecular pump and a rotary vane pump with a filling of special oil may be supplied
- For especially sensitive processes also a dry compressing vacuum pump like the ECODRY M may be used as the backing pump

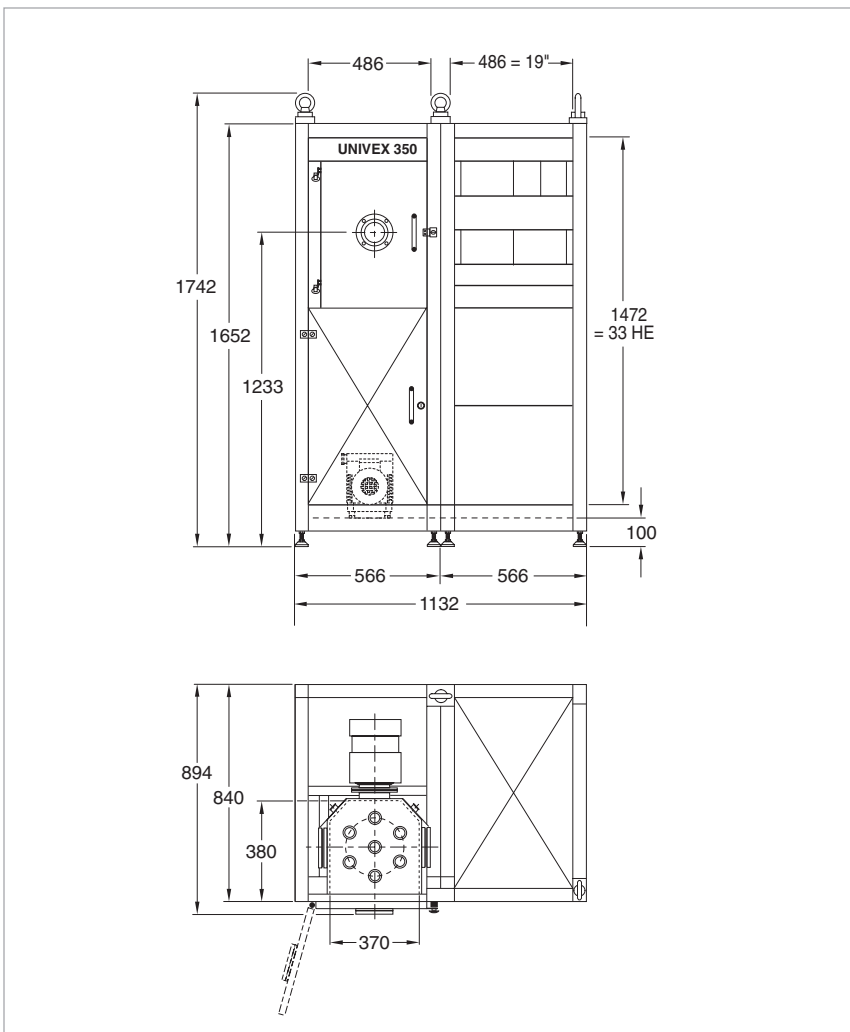

Dimensional drawing for the UNIVEX 350

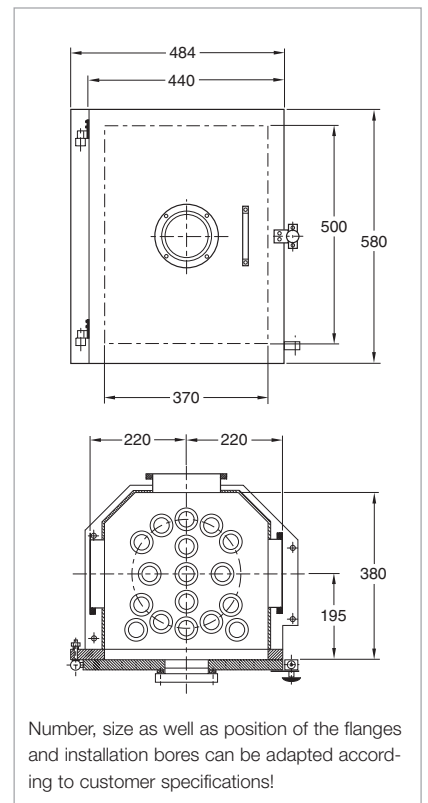

Number, size as well as position of the flanges and installation bores can be adapted according to customer specifications!

Dimensional drawing for the vacuum chamber

Technical Data

UNIVEX 350

Vacuum chamber		
Material		Stainless steel
Dimensions		
Inside width	mm	370
Inside depth	mm	380
Inside height	mm	500
Connections ¹⁾		
Front side		Door with window
Rear side	DN	160 ISO-K (pump system connection), 2 x 16 ISO-KF, 2 x 40 ISO-KF
Bottom plate		15 installation holes 34.5 mm dia.
Cover plate		7 installation holes 34.5 mm dia.
Left side	DN	160 ISO-K
Right side	DN	160 ISO-K
High vacuum pump		TURBOVAC 600 C
Nominal pumping speed for N ₂	l x s ⁻¹	560
Power supply		TURBOTRONIK NT 20
Backing pump		TRIVAC D 25 B
Nominal pumping speed	m ³ x h ⁻¹	29.5
Controller		PLC with graphic touchscreen
Required supplies		
Voltage		400 V, 3 phases + N, 50/60 Hz ²⁾
Cooling water		
Inlet pressure	bar (abs.)	4 to 7
Consumption, approx.	l x min ⁻¹	1 ³⁾
Feed temperature	°C	15 to 25
Compressed air	bar (abs.)	4 to 7
Weight, approx.	kg	400 ³⁾

Ordering Information

UNIVEX 350

UNIVEX 350	upon request
------------	--------------

¹⁾ Standard configuration, other hole patterns/flanges upon request

²⁾ Other voltages and frequencies upon request

³⁾ Without chamber installations/process components

UNIVEX 350 G



UNIVEX 350 G, consisting of electrical cabinet (left) and coating module (right)

Laboratory System for Glove Box Attachment with Cubic Vacuum Chamber and Front Sliding Door (350 mm diameter)

Typical Applications

- Vacuum coating in research and development
- Pre-production trials
- Special experiments

Design

- The UNIVEX 350 G consists of a coating module and a 19" rack cabinet
- Installed in the coating module is the process chamber and the pump system
- Accommodated in the cabinet is the electric power supply with the PLC based system controller with graphic touchscreen. Moreover, the supply units for the process components are also accommodated in this cabinet section

Vacuum Chamber

- The base plate is attached to the base frame
- The sliding front door is operated manually and closed pneumatically
- Bottom plate and lid are provided with installation holes
- Additional flanges at the side for pump system and process components

Pump System

- The standard pumping equipment comprises a TRIVAC D 25 B two-stage rotary vane pump and a TURBOVAC 600 C turbomolecular pump

Vacuum Measurement


- Depending on the type of application, a combination vacuum gauge operating according to the cold cathode or hot cathode principle may be installed

Advantages to the User

- Intended for fitting to the rear side of a glove box
- Convenient access to the process through the glove box by means of the sliding front door
- Easy access to the chamber unit through the rear service door
- Only the sliding door is located in the glove box. All other system components can be easily accessed from the outside
- Modular system design
- Freedom of installation and retrofitting of process components without problems
- Very simple to operate and use via programmable control
- For RF sputtering
- Pump system adapted to the individual process

Options

- Easy to disassemble coating protection panels are available
- For short pumpdown times, a bypass line bypassing the high vacuum pump can be provided
- For processes producing increased amounts of gas or for low operating pressures, the UNIVEX 350 G may also be equipped with turbomolecular pumps offering a higher pumping speed (TURBOVAC 1100 C, for example) or cryogenic pumps may be specified
- For processes in which aggressive media need to be pumped, a seal gas version of the turbomolecular pump and a rotary vane pump with a special oil filling can be supplied
- For particularly sensitive processes, a dry compressing backing pump like the SCROLLVAC SC 30 D can be used

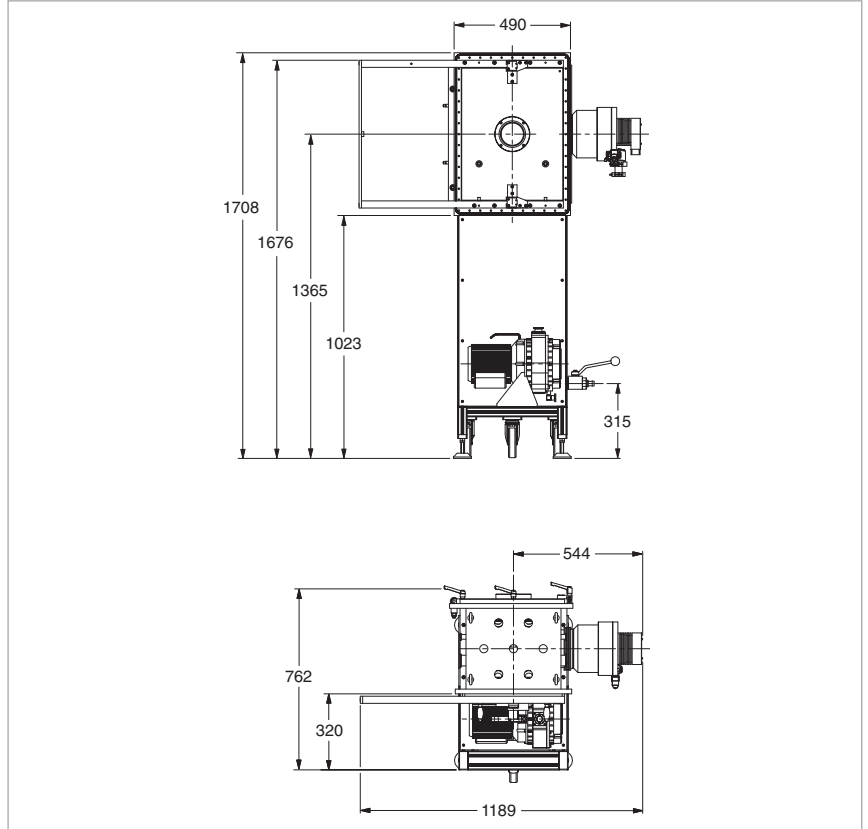


Dimensional drawing for the UNIVEX 350 G

Turnkey Solutions

Upon request Oerlikon Leybold Vacuum will also arrange the delivery of turnkey solutions consisting of the UNIVEX 350 G coating system and a glove box from a single source.

More information on such systems is available upon request.

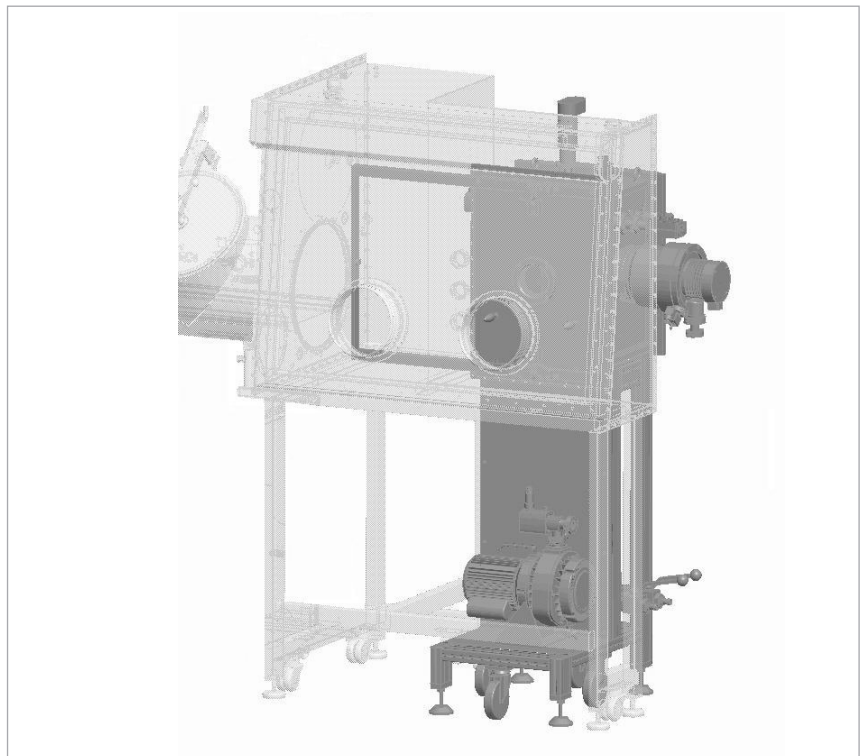

3D view of a glove box with the UNIVEX 350 G coating module fitted to the rear

Technical Data

UNIVEX 350 G

Vacuum chamber		
Material		Stainless steel
Dimensions		
Inside width	mm	370
Inside depth	mm	380
Inside height	mm	500
Connections ¹⁾		
Front side		Sliding door for glove box access; manually operated and pneumatically closing
rear side		turning door for service access; manually locked
Bottom plate		15 installation holes 34.5 mm dia.
Cover plate		7 installation holes 34.5 mm dia.
Left side	DN	160 ISO-K (pump system connection), 40 ISO-KF, 16 ISO-KF
Right side	DN	40 ISO-KF, 16 ISO-KF
High vacuum pump		TURBOVAC 600 C
Nominal pumping speed for N ₂	l x s ⁻¹	560
Power supply		TURBOTRONIK NT 20
Backing pump		TRIVAC D 25 B
Nominal pumping speed	m ³ x h ⁻¹	29.5
Controller		PLC with graphic touchscreen
Required supplies		
Voltage		400 V, 3 phases + N, 50/60 Hz ²⁾
Cooling water		
Inlet pressure	bar (abs.)	4 to 7
Consumption, approx.	l x min ⁻¹	1 ³⁾
Feed temperature		15 to 25 °C
Compressed air	bar (abs.)	4 to 7
Weight, approx.	kg	350 ³⁾

Ordering Information

UNIVEX 350 G

UNIVEX 350 G	upon request
--------------	---------------------

¹⁾ Standard configuration, other hole patterns/flanges/viewing windows upon request


²⁾ Other voltages and frequencies upon request

³⁾ Applies only to the coating module, without chamber units/process components

UNIVEX 450 B



UNIVEX 450 B with water-cooled vacuum chamber (option)

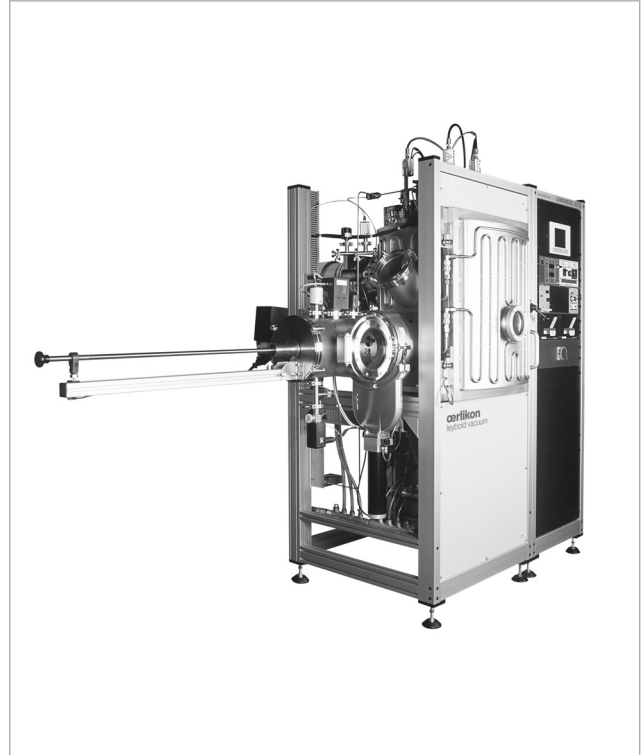

UNIVEX 450 B with water-cooled process chamber and manually operated load lock chamber (option)

Laboratory System with Cubic Vacuum Chamber and Front Door (500 mm diameter)

Typical Applications

- Vacuum coating in research and development
- Pre-production trials
- Special experiments

Design

- The UNIVEX 450 B consists of two separable 19" cabinet sections
- Accommodated in one cabinet section is the process chamber and the pump system
- Accommodated in the second cabinet section is the electric power supply with the PLC based system controller with graphic touchscreen. Moreover, the supply units for the process components are also accommodated in this cabinet section

Vacuum Chamber

- The base plate is attached to the base frame
- The door is equipped with a viewing window
- Bottom plate and lid are provided with installation holes
- Additional flanges at the side for pump system and process components

Pump System

- The standard pumping equipment comprises a TRIVAC D 65 B two-stage rotary vane pump and a TURBOVAC 1100 C turbomolecular pump

Vacuum Measurement


- Depending on the type of application, a combination vacuum gauge operating according to the cold cathode or hot cathode principle may be installed

Advantages to the User


- Modular system design
- Freedom of installation and retrofitting of process components without problems
- Vacuum chamber with door
- Convenient access to the chamber installations
- Very simple to operate and use via programmable control
- Suited for being installed in a clean room wall
- For RF sputtering
- Pump system adapted to the individual process

Options

- The chamber can also be delivered in a water-cooled version
- Easy to disassemble coating protection panels are available
- For short pumpdown times, a bypass line bypassing the high vacuum pump can be provided
- Fitting of a second coating module with a vacuum chamber (fitted to the right of the electrical cabinet) is provided for
- For processes producing increased amounts of gas or for low operating pressures, the UNIVEX 450 B may also be equipped with turbomolecular pumps offering a higher pumping speed (TURBOVAC T1600, for example) or cryogenic pumps may be specified
- For processes in which aggressive media need to be pumped, a seal gas version of the turbomolecular pump and a rotary vane pump with a special oil filling can be supplied
- For particularly sensitive processes, a dry compressing backing pump like the SCROLLVAC SC 30 D can be used

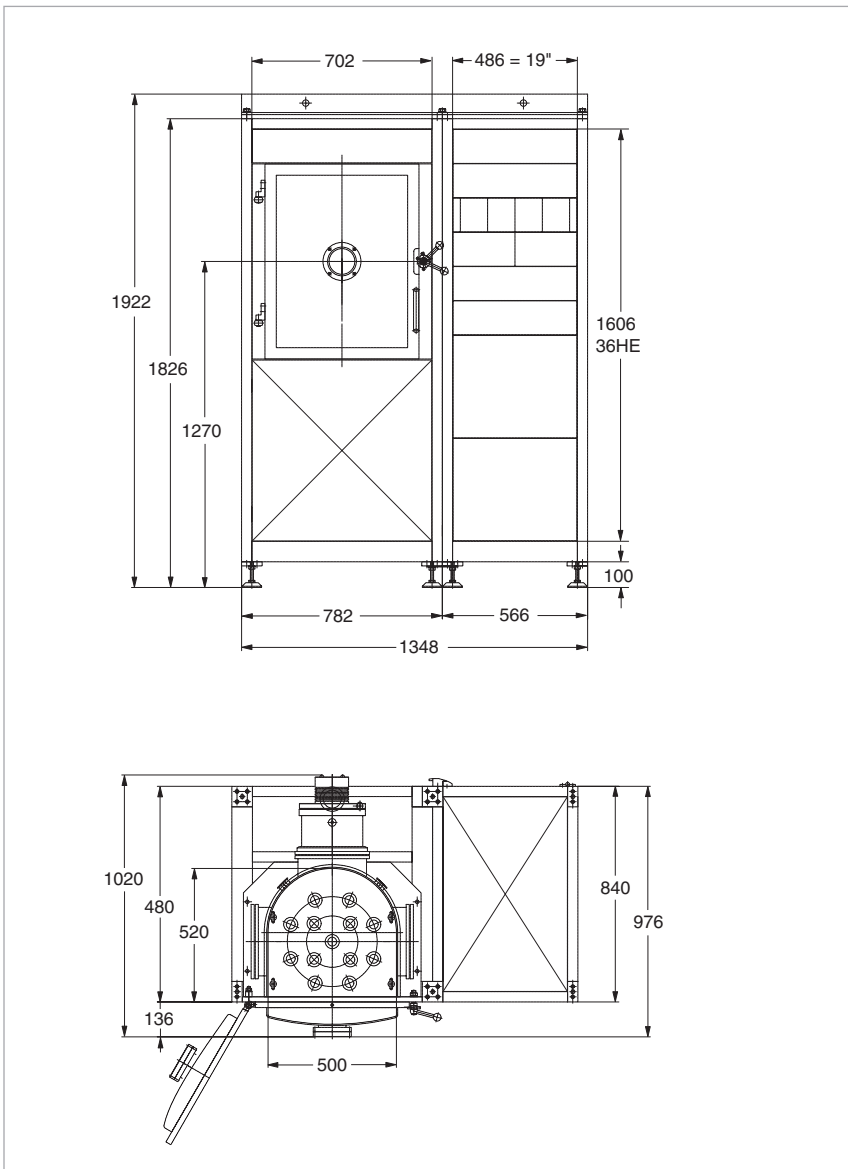

Dimensional drawing for the UNIVEX 450 B

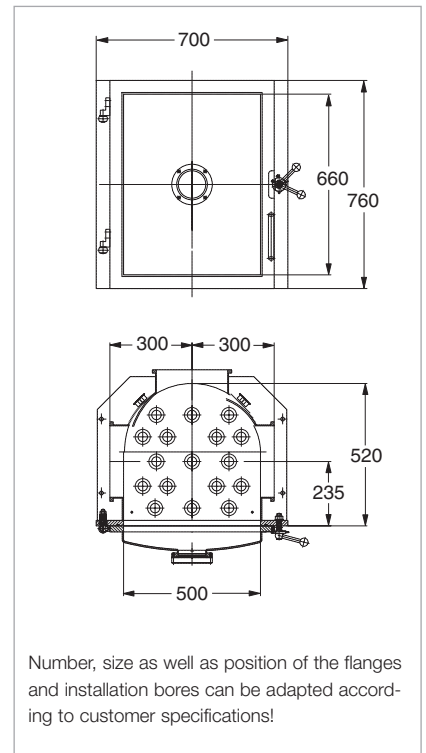

Number, size as well as position of the flanges and installation bores can be adapted according to customer specifications!

Dimensional drawing for the vacuum chamber

Technical Data

UNIVEX 450 B

Vacuum chamber		
Material		Stainless steel
Dimensions		
Inside width	mm	500
Inside depth	mm	500
Inside height	mm	650
Connections ¹⁾		
Front side		Door with window
Rear side	DN	250 ISO-K (pump system connection), 4 x 16 ISO-KF, 2 x 40 ISO-KF
Bottom plate		20 installation holes 34.5 mm dia.
Cover plate		10 installation holes 34.5 mm dia.
Left side	DN	250 ISO-K
Right side	DN	250 ISO-K
High vacuum pump		TURBOVAC 1100 C
Nominal pumping speed for N ₂	l x s ⁻¹	1050
Power supply		TURBOTRONIK NT 20
Backing pump		TRIVAC D 25 B
Nominal pumping speed	m ³ x h ⁻¹	29.5
Controller		PLC with graphic touchscreen
Required supplies		
Voltage		400 V, 3 phases + N, 50/60 Hz ²⁾
Cooling water		
Inlet pressure	bar (abs.)	4 to 7
Consumption, approx.	l x min ⁻¹	1 ³⁾
Feed temperature	°C	15 to 25
Compressed air	bar (abs.)	4 to 7
Weight, approx.	kg	500 ³⁾

Ordering Information

UNIVEX 450 B

UNIVEX 450 B

upon request

¹⁾ Standard configuration, other flanges/hole patterns/viewing windows upon request

²⁾ Other voltages and frequencies upon request

³⁾ Applies only to the basic unit without coating equipment

UNIVEX 450 for Dactyloscopy



UNIVEX 450 for dactyloscopy

Dactyloscopy is a term from the area of criminal investigation meaning: "Identification of a person through his fingerprints". Depending on the material of the part which was touched, different methods are used to render the fingerprints visible.

In the case of materials like plastic shopping bags, for example, foils, handlebars etc. evaporation methods have been found to be most useful.

The method itself utilises the effect well known from normal evaporation processes where the evaporated material will adhere better (and thicker) on the skin material (water, amino acids, fat and alike) deposited by the finger compared to the surrounding untouched material.

An optimum contrast is attained by selecting a suitable evaporation material, usually gold or zinc.

Benefits of this method

- No "smearing" of existing traces compared to conventional methods
- Large surface areas (up to 80 x 40 cm max.) carrying fingerprints can be checked completely in one pass
- The time needed for one pass is only about 10 minutes (depending on the material carrying the fingerprints)
- Good contrast also in the case of multicolour surfaces
- Fixation of the deposited material with the traces is easy – the results may be well documented (can be photographed)
- The carrier of the fingerprints is not destroyed

UNIVEX 450 C



For special applications we can also supply cluster systems based on the UNIVEX concept. These clusters are equipped according to customers requirements and incorporate separate processing and load lock and transfer chambers.

UNIVEX 450 C with coating module and electrical cabinet (example photograph).
The coating module consists of two process chambers (left and right) as well as the loadlock and transfer chamber in between

Test Systems with a Vacuum Chamber



We can also supply vacuum chambers with custom pump systems for testing of various components.

Test system with a 700 mm dia chamber

Accessories

Standard Accessories for UNIVEX Systems

Blank-Off Screw Fitting

For 34.5 mm dia. hole.



Blank-off screw fitting

PS 113 A Safety Switch

For safety interlock arrangements in connection with the UNIVEX system controller, respectively optionally connected power supply equipment (for sputtering, electron beam evaporation or vacuum etching, for example).



PS 113 A Safety Switch

Variable Leak Valve with Isolation Valve

For manually controlled admission of gas in connection with plasma processes (sputtering, etching and glow discharge cleaning).

Technical Data

Material	Stainless steel	
Seal	FPM	
Weight	kg	0.1

Ordering Information

Blank-off screw fitting	Part No. 030 40
-------------------------	------------------------

Blank-Off Screw Fitting

Blank-Off Screw Fitting

Technical Data

Switching pressure	mbar	approx. 6 below atmospheric pressure
Return switching pressure	mbar	3 below atmospheric pressure
Switching inaccuracy	mbar	26
Max. permissible operating pressure (abs.)	mbar	2000
Storage temperature range	°C	-25 to +85
Nominal temperature range	°C	0 to 85
Switching contact	Changeover contacts, gold-plated, for prog. controls	
Contact life	> 10 ⁵ switching cycles	
Switching capacity	100 mA / 24 V AC 30 mA / 24 V DC	
Electrical connection	6.3 mm flat plug	
Vacuum connection	DN	16 ISO-KF
Materials in contact with the medium	Stainless steel 1.4305, 1.4310, Stainless steel 1.4300 PTFE coated	
Protection class	IP	44

Ordering Information

Low pressure safety switch
PS 113 A, DN 16 ISO-KF;
complete with 3 m long cable

Safety Switch

Part No. 230 011

Technical Data

Gas admission rate q_L	mbar x l x s ⁻¹	5 x 10 ⁻⁶ to 1 x 10 ³
Connection flange	DN	16 ISO-KF

Ordering Information

Variable leak valve with isolation valve	Part No. 215 010
--	-------------------------

Variable Leak Valve with Isolation Valve

Variable Leak Valve with Isolation Valve

(see also Product Section C14 "Vacuum Valves")

Gas Flow Controller (MFC)


For controlled admission of gas in connection with automated plasma processes (sputtering, etching, glow discharge cleaning). Intended to be remotely controlled by the customer's PC or PLC, i.e. a separate MFC controller unit is not included in the delivery.

Electromotively Operated Vapor Source Shutter

For covering the source during thermal or electron beam evaporation. With gear motor and shutter panel; can be fitted to the inside sections of the chamber.

6-Way Measurement Feedthroughs

For connection of the vapor source shutter; for 34.5 mm holes, plug-in soldered contact on the inside.

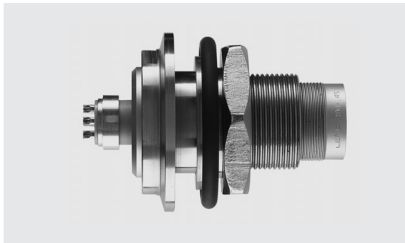

6-way measurement feedthrough

Control Cable, 6-Way

For connection between measurement feedthrough and power supply unit for the vapor source shutter, complete with connection plugs.

Interlocking Kit

For providing touch protection against high-voltage carrying parts within the vacuum chamber. The safety contact must be connected to the interlock input at the high-voltage power supply unit thereby ensuring that the equipment can be enabled only while the chamber door is closed.

Technical Data

Gas flow, max.	sccm	selectable between 10 and 500
Supply voltage	V DC	24
Control signal	V	0 - 5 analog

Ordering Information

Gas flow controller	upon request
---------------------	---------------------

Gas Flow controller

Technical Data

Control voltage	V DC	24 (pulses per second)
Dimensions of the shutter panel	mm	different, for example 42 x 42 or 100 mm dia.
Weight	kg	0.2

Ordering Information

Electromotively operated vapor source shutter	upon request
---	---------------------

Vapor Source Shutter

Technical Data

Rating per conductor	V A	max. 700 16
Seal		FPM
Weight	kg	0.3

Ordering Information

6-way measurement feedthrough	upon request
-------------------------------	---------------------

Measurement Feedthrough

Technical Data

Length	m	3
Weight	kg	0.2

Ordering Information

6-way control cable	upon request
---------------------	---------------------

Control Cable

Technical Data

Mechanical closing contact	Safety door switch
Electrical closing contact	floating

Ordering Information

Interlocking kit for UNIVEX vacuum chamber	upon request
--	---------------------

Interlocking Kit

Interlocking Kit

Components for Thermal Evaporation of High Melting Point Materials (metals)

Single Thermal Evaporator

Consisting of two water-cooled high voltage feedthroughs with terminal blocks for 34.5 mm dia. holes.



Single thermal evaporator

Dual Thermal Evaporator

Consisting of three water-cooled high voltage feedthroughs with terminal blocks for 34.5 mm dia. holes.



Dual thermal evaporator

Power Supply Cables

For single and dual thermal evaporators, equipped with terminals and clamping pieces.

Technical Data

Rating per conductor	V	max. 100
	A	500
Seals		FPM
Water connection	mm	Hose 4/6 dia.2.5
Weight	kg	

Single Thermal Evaporator

Ordering Information

Single thermal evaporator

Single Thermal Evaporator

upon request

Technical Data

Rating per conductor	V	max. 100
	A	500
Seals		FPM
Water connection	mm	Hose 4/6 dia.
Weight	kg	3.9

Dual Thermal Evaporator

Ordering Information

Dual thermal evaporator

Dual Thermal Evaporator

upon request

Technical Data

Length	m	2 ¹⁾
Rating	V	max. 100
	A	500
Cross section	mm ²	120
Weight	kg	3.5

Power Supply Cable

Ordering Information

Power supply cable

Power Supply Cable

upon request ²⁾


¹⁾ Standard length. Other lengths can be specified

²⁾ For the single thermal evaporator, two supply lines are required.

For the dual thermal evaporator, three supply lines are required

AS 053 Power Supply Unit

For supplying thermal evaporators and one solenoid-actuated source shutter. With LCD display for current read out; with membrane key pad.

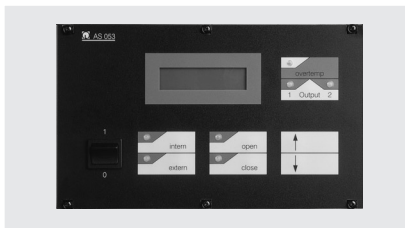

AS 053 power supply unit

AS 053/2 Power Supply Unit

For supplying power to two thermal evaporators with vapor source shutters. With LCD display for current read out; with membrane key pad.

Technical Data

Cabinet	mm	1/2 19" rack module, 3 HU 400 deep
Outputs		1 x evaporator output, 5 V, 400 A max. can be rewired to 10 V, 200 A max. 1 x shutter output, 24 V DC, 1 s pulse
Inputs		Remote control unit for controlling the evaporation power (0 to 10 V) Remote control for the shutter
Main power supply		230 V, 50/60 Hz, 10 A
Weight	kg	15

Power Supply Unit

Ordering Information

AS 053 power supply unit	upon request
--------------------------	---------------------

Power Supply Unit

Technical Data

Cabinet	mm	19" rack module, 3 HU 400 deep
Outputs		2 x evaporator output, 5 V, 400 A max. can be rewired to 10 V, 200 A max. 2 x shutter output, 24 V DC, 1 s pulse
Inputs		Remote control unit for controlling the evaporation power (0 to 10 V) Remote control for the shutter Switchover evaporator 1 / 2
Main power supply		230 V, 50/60 Hz, 10 A
Weight	kg	30

Power Supply Unit

Ordering Information

AS 053/2 power supply unit	upon request
----------------------------	---------------------

Power Supply Unit

Components for Thermal Evaporation of Low Melting Point (organic) Materials

For the purpose of evaporating temperature sensitive materials, commonly of an organic nature, Oerlikon Leybold Vacuum is offering special organic material evaporators.

These ensure a coating process at precisely controlled evaporation temperatures which typically range between 200 °C and 400 °C.

For installation within the UNIVEX systems, Oerlikon Leybold Vacuum supplies organic material evaporators as a complete package, consisting of evaporator source, automatic shutter, rotary vacuum feed through and 19" rack mount controller.

Components for Electron-Beam Evaporation

Various models of electron-beam evaporators and power supplies are available for installation in the UNIVEX systems.

Electron-Beam Evaporator

The selection of a suitable electron beam evaporator will chiefly depend on the amount of available space, the desired evaporation rate and the film thickness as well as the number and type of materials which need to be evaporated. Single crucible as well as rotatable multi-crucible evaporators are available.

Power Supplies

The power supply unit for the individual electron beam evaporators is selected depending on the maximum evaporation power which is required, as well as the demanded properties for X/Y beam deflection. Models with output power ratings ranging from 3 kW to 10 kW are available.

As a rule, the maximum output power of the power supply may not exceed the maximum permissible power specified for the evaporator.

Safety Regulations for Electron-Beam Evaporator Applications

When installing electron beam evaporators within the UNIVEX 300, only the stainless steel bell jar must be used. Moreover, a safety interlocking system will be necessary for all UNIVEX types. For the UNIVEX 300 a separate interlocking kit is available; in the case of the UNIVEX 350 and 450 B this kit is already included.

As further safety means a water flow monitor is required for each electron-beam evaporation unit so as to ensure intensive cooling of the electron-beam evaporator.

This water flow monitor is included in the delivery of a UNIVEX system equipped with an electron beam evaporator.

Upon request we shall be pleased to provide an offer which specifically matches the requirements of your application.

Components for Sputtering

DC Sputtering

Widely differing DC sputtering sources may be fitted within all UNIVEX systems. Their selection will depend on the size of the substrate, the required target material and the available installation space. Circular planar sputtering sources of 50 mm to 200 mm in diameter as well as various rectangular sources are available. The power supply units (providing an output power between 500 W and 3 kW) may be installed within the UNIVEX 19" electrical cabinets.

DC sputtering equipment is suited for all UNIVEX systems.

RF Sputtering

Widely differing RF sputtering sources may be fitted within all UNIVEX systems (exception: UNIVEX 300). Their selection will depend on the size of the substrate, the required target material and the available installation space. Circular planar sputtering sources of 50 mm to 200 mm in diameter as well as various rectangular sources together with the necessary RF matching components are available. The power supply units (providing an output power between 150 W and 1.5 kW) may be installed within the UNIVEX systems.

Gas Inlet

Sputtering sources can only be operated with gas admission. For this, manually operated variable leak valves up to automatically controlled mass flow controllers are available.

Throttling the Pumping Speed

In order to protect the high vacuum pump against the high process pressure present during plasma operation, and to reduce process gas consumption, the UNIVEX systems are generally equipped with high vacuum gate valves having three positions which are fitted between chamber and high vacuum pump.

Safety Regulations in Connection with Sputtering Applications

When installing electron beam evaporators within the UNIVEX 300, only the stainless steel bell jar must be used. Moreover, a safety interlocking system will be necessary for all UNIVEX types. For the UNIVEX 300 a separate interlocking kit is available; in the case of the UNIVEX 350 and 450 B this kit is already included.

As further safety means a water flow monitor is required for each electron-beam evaporation unit so as to ensure intensive cooling of the electron-beam evaporator.

This water flow monitor is included in the delivery of a UNIVEX system equipped with an electron beam evaporator.

Upon request we shall be pleased to provide an offer which specifically matches the requirements of your application.

Components for Glow Discharge Cleaning

With glow discharge electrode, high voltage feedthrough for 34.5 mm dia. hole and connection cable for fitting to the central rotary feedthrough.



Glow discharge assembly with high voltage feedthrough and connection cable

Technical Data

Electrode material		Aluminum
Insulation		Ceramics
Max. ratings	V mA	2000 65
Sealing material of the high voltage feedthrough		FPM
Length of the connection cable	m	2
Weight	kg	1

Glow Discharge Assembly

C 2000 High Voltage Power Supply Unit

For supplying the glow discharge assembly.



C 2000 high voltage power supply unit

Ordering Information

Glow discharge assembly for UNIVEX 300, 350 and 350 G for UNIVEX 450 B	upon request upon request
--	--

Glow Discharge Assembly

Technical Data

Cabinet		19" rack module, 3 HU
Output	V mA	max. 2000 max. 65
Electrical connection		230 V, 50/60 Hz, 150 VA
Remote control and locking input		included

High Voltage Power Supply Unit

Ordering Information

C 2000 high voltage power supply unit	upon request
---------------------------------------	---------------------

High Voltage Power Supply Unit

Gas Admission

Operation of the glow discharge cleaning facility is only possible by admitting gas. For this, manually operated dosing valves up to automatically controlled mass flow controllers offer all possibilities.

Throttling the Pumping Speed

In order to protect the high vacuum pump against the high process pressure present during plasma operation, and reduce process gas consumption, the UNIVEX systems are generally equipped with high vacuum gate valves having three positions which are fitted between chamber and high vacuum pump.

Components for Film Thickness Measurements


Various thin film thickness measuring instruments may be installed in the UNIVEX units.

The selection depends on the demanded measurements tasks and the required degree of automation.

As standard, oscillating crystal systems are used. These may consist of one or several sensing heads with or without shutter, and upon request are available also for UHV operation (i.e. are suitable for degassing).

These are driven either by a monitor (allowing only the measurement of deposition rate and film thickness) or by a controller (allowing besides measurement of the film parameters also to control the deposition rate).

Upon request we shall be pleased to provide an offer which specifically matches the requirements of your application.

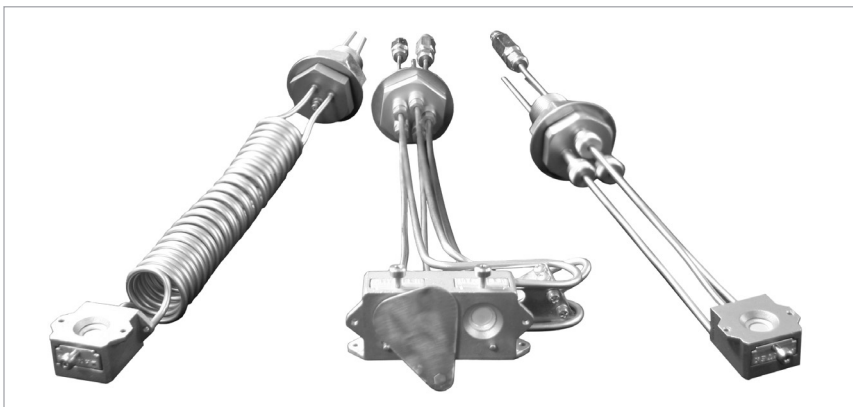


Substrate Manipulation Solutions

Substrate Rotation

In order to improve the film properties in connection with deposition processes in a vacuum, a rotary movement of the substrates will often be necessary.

For this reason Oerlikon Leybold Vacuum is offering for the UNIVEX product family a broad spectrum of possible solutions. These begin with a simple manually operated rotatable axle and end with speed controlled planetary gears. Moreover, there exists the possibility of linking the rotary movement to additional features like heating, cooling or bias. Depending on the number of substrates and their dimensions as well as the temperature requirements, a check of the specific circumstances will be necessary in each case.


More information on this is available upon request.

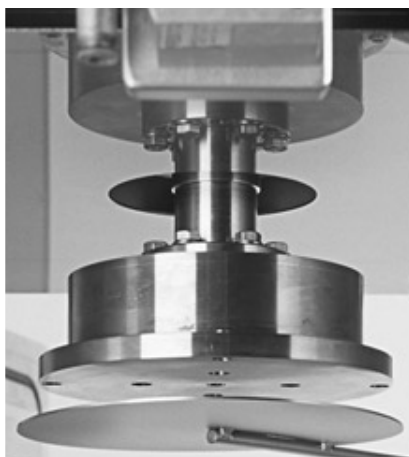

Rotatable, temperature controlled substrate holder with substrate shutter

Substrate Holding

For mounting the substrate within the chamber, Oerlikon Leybold Vacuum offers vacuum substrate holders manufactured according to customer's specifications.

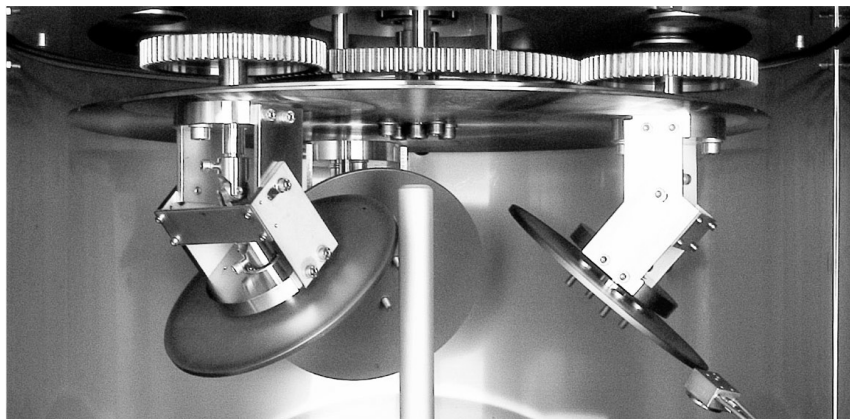

Planetary gear

Substrate Heating

For temperature controlled heating of substrates, Oerlikon Leybold Vacuum offers different processes (resistance heaters, quartz lamp heaters etc., for example). Depending on the size of the substrate and the temperature range, we are prepared to offer a customised heating solution.

Substrate Cooling

Sensitive substrates or substrate masks require intensive cooling while the coating process is in progress. For this, Oerlikon Leybold Vacuum delivers substrate holders for different temperature ranges (cold water, liquid nitrogen or special fluids).

Substrate Bias

For the purpose of cleaning substrates before coating them or for improving the adhesive properties during the coating process, RF or DC biasing of the substrates can provide the required advantage. For this, Oerlikon Leybold Vacuum can supply insulated substrate mounts with or without the corresponding supply units.

Sales and Service

Germany

Oerlikon
Leybold Vacuum GmbH
Bonner Strasse 498
D-50968 Cologne
Phone: +49-(0)221-347 1234
Fax: +49-(0)221-347 1245
sales.vacuum@oerlikon.com
www.oerlikon.com

Oerlikon
Leybold Vacuum GmbH
Sales Area North/Northeast
Branch Office Berlin
Buschkrugallee 33
1. Obergeschoss
D-12359 Berlin
Phone: +49-(0)30-435 609 0
Fax: +49-(0)30-435 609 10
sales.vacuum.bn@oerlikon.com

Oerlikon
Leybold Vacuum GmbH
Sales Area South/Southwest
Branch Office Munich
Sendlinger Strasse 7
D-80331 Munich
Phone: +49-(0)89-357 33 9-10
Fax: +49-(0)89-357 33 9-33
sales.vacuum.mn@oerlikon.com
service.vacuum.mn@oerlikon.com

Oerlikon
Leybold Vacuum GmbH
Sales Area West & Benelux
Branch Office Cologne
Bonner Strasse 498
D-50968 Cologne
Phone: +49-(0)221-347 1270
Fax: +49-(0)221-347 1291
sales.vacuum.kn@oerlikon.com

Oerlikon
Leybold Vacuum GmbH
Service Competence Center
Emil-Hoffmann-Strasse 43
D-50996 Cologne-Suerth
Phone: +49-(0)221-347 1439
Fax: +49-(0)221-347 1945
service.vacuum.kn@oerlikon.com

Oerlikon
Leybold Vacuum GmbH
Mobil Customer Service
Emil-Hoffmann-Strasse 43
D-50996 Cologne-Suerth
Phone: +49-(0)221-347 1765
Fax: +49-(0)221-347 1944
service.vacuum.kn@oerlikon.com

Oerlikon
Leybold Vacuum GmbH,
Dresden
Zur Wetterwarte 50, Haus 304
D-01109 Dresden
Service:
Phone: +49-(0)351-88 55 00
Fax: +49-(0)351-88 55 041
info.vacuum.dr@oerlikon.com

Oerlikon
Leybold Vacuum USA Inc.
5700 Mellon Road
USA-Export, PA 15632
Phone: +1-724-327-5700
Fax: +1-724-325-3577
info.vacuum.ex@oerlikon.com

Europe

Belgium
Oerlikon
Leybold Vacuum Nederland B.V.
Belgisch bijkantoor
Leuvensesteenweg 542-9A
B-1930 Zaventem
Sales:
Phone: +32-2-711 00 83
Fax: +32-2-720 83 38
sales.vacuum.zv@oerlikon.com
Service:
Phone: +32-2-711 00 82
Fax: +32-2-720 83 38
service.vacuum.zv@oerlikon.com

France
Oerlikon
Leybold Vacuum France S.A.
7, Avenue du Québec
Z.A. Courtaboef 1 - B.P. 42
F-91942 Courtaboef Cedex
Sales and Service:
Phone: +33-1-69 82 48 00
Fax: +33-1-69 07 57 38
sales.vacuum.or@oerlikon.com

Oerlikon
Leybold Vacuum France S.A.
Valence Factory
640, Rue A. Bergès - B.P. 107
F-26501 Bourg-lès-Valence Cedex
Phone: +33-4-75 82 33 00
Fax: +33-4-75 82 92 69
info.vacuum.vc@oerlikon.com

Great Britain
Oerlikon
Leybold Vacuum UK Ltd.
Unit 2
Silverglade Business Park
Leatherhead Road
UK-Chessington, Surrey KT9 2QL
Sales:
Phone: +44-13-7273 7300
Fax: +44-13-7273 7301
sales.vacuum.ln@oerlikon.com
Service:
Phone: +44-20-8971 7030
Fax: +44-20-8971 7003
service.vacuum.ln@oerlikon.com

Italy
Oerlikon
Leybold Vacuum Italia S.p.A.
8, Via Trasimeno
I-20128 Milano
Sales:
Phone: +39-02-27 22 31
Fax: +39-02-27 20 96 41
sales.vacuum.mi@oerlikon.com
Service:
Phone: +39-02-27 22 31
Fax: +39-02-27 22 32 17
service.vacuum.mi@oerlikon.com

Oerlikon
Leybold Vacuum Italia S.p.A.
Field Service Base
Z.I. Le Capanne
I-05021 Acquasparta (TR)
Phone: +39-0744-93 03 93
Fax: +39-0744-94 42 87
service.vacuum.mi@oerlikon.com

Netherlands
Oerlikon
Leybold Vacuum Nederland B.V.
Proostwetering 24N
NL-3543 AE Utrecht
Sales and Service:
Phone: +31-(30) 242 6330
Fax: +31-(30) 242 6331
sales.vacuum.ut@oerlikon.com
service.vacuum.ut@oerlikon.com

Spain
Oerlikon
Leybold Vacuum Spain, S.A.
C/. Huelva, 7
E-08940 Cornellà de Llobregat
(Barcelona)
Sales:
Phone: +34-93-666 46 16
Fax: +34-93-666 43 70
sales.vacuum.ba@oerlikon.com
Service:
Phone: +34-93-666 49 51
Fax: +34-93-685 40 10
service.vacuum.ba@oerlikon.com

Sweden
Oerlikon
Leybold Vacuum Scandinavia AB
Box 9084
SE-40092 Göteborg
Sales and Service:
Phone: +46-31-68 84 70
Fax: +46-31-68 39 39
info.vacuum.gt@oerlikon.com
Visiting/delivery address:
Datavägen 57B
SE-43632 Askim

Switzerland
Oerlikon
Leybold Vacuum Schweiz AG
Leutschenbachstrasse 55
CH-8050 Zürich
Sales:
Phone: +41-044-308 40 50
Fax: +41-044-302 43 73
sales.vacuum.zh@oerlikon.com
Service:
Phone: +41-044-308 40 62
Fax: +41-044-308 40 60
service.vacuum.zh@oerlikon.com

America

USA
Oerlikon
Leybold Vacuum USA Inc.
5700 Mellon Road
USA-Export, PA 15632
Phone: +1-724-327-5700
Fax: +1-724-325-3577
info.vacuum.ex@oerlikon.com
Sales:
Eastern & Central time zones
Phone: +1-724-327-5700
Fax: +1-724-333-1217
Pacific, Mountain, Alaskan &
Hawaiian time zones
Phone: +1-480-752-9191
Fax: +1-480-752-9494
Service:
Phone: +1-724-327-5700
Fax: +1-724-325-3577

Oerlikon
Leybold Vacuum GmbH
Bonner Strasse 498
D-50968 Cologne
Phone: +49-(0)221-347 0
Fax: +49-(0)221-347 1250
info.vacuum@oerlikon.com

Asia

P.R. China
Oerlikon
Leybold Vacuum (Tianjin)
International Trade Co. Ltd.
Beichen Economic
Development Area (BEDA),
Shanghai Road
Tianjin 300400
China
Sales and Service:
Phone: +86-22-2697 0808
Fax: +86-22-2697 4061
Fax: +86-22-2697 2017
sales.vacuum.tj@oerlikon.com
service.vacuum.tj@oerlikon.com

Oerlikon
Leybold Vacuum
(Tianjin) Co. Ltd.
Beichen Economic
Development Area (BEDA),
Shanghai Road
Tianjin 300400
China
Sales and Service:
Phone: +86-22-2697 0808
Fax: +86-22-2697 4061
Fax: +86-22-2697 2017
info.vacuum.tj@oerlikon.com
service.vacuum.tj@oerlikon.com

Oerlikon
Leybold Vacuum (Tianjin)
International Trade Co. Ltd.
Shanghai Branch:
Add: No.33
76 Futedong San Rd.
Waigaoqiao FTZ
Shanghai 200131
China
Sales and Service:
Phone: +86-21-5064-4666
Fax: +86-21-5064-4668
info.vacuum.sh@oerlikon.com
service.vacuum.tj@oerlikon.com

Oerlikon
Leybold Vacuum (Tianjin)
International Trade Co. Ltd.
Guangzhou Office and
Service Center
1st F, Main Building
Science City Plaza,
No.111 Science Revenue,
Guangzhou Science City
(GZSC) 510663, Guangzhou,
China
Sales:
Phone: +86-20-8723-7873
Phone: +86-20-8723-7597
Fax: +86-20-8723-7875
info.vacuum.gz@oerlikon.com
service.vacuum.tj@oerlikon.com

Oerlikon
Leybold Vacuum (Tianjin)
International Trade Co. Ltd.
Beijing Branch:
1-908, Beijing Landmark Towers
8 North Dongsanhuan Road
Chaoyang District
Beijing 100004
China
Sales:
Phone: +86-10-6590-7622
Fax: +86-10-6590-7607
sales.vacuum.bj@oerlikon.com

India
Oerlikon
Leybold Vacuum India Pvt Ltd.
EL-22, J Block
MIDC Bhosari
Pune 411026
India
Sales and Service:
Phone: +91-20-3061 60000
Fax: +91-20-2712 1571
sales.vacuum.pu@oerlikon.com
service.vacuum.pu@oerlikon.com

Japan
Oerlikon
Leybold Vacuum
Japan Co., Ltd.
Headquarter
23-3, Shin-Yokohama
3-chome
Tobu A.K. Bldg. 4th Floor
Kohoku-ku
Yokohama-shi 222-0033
Sales:
Phone: +81-45-471-3330
Fax: +81-45-471-3323
info.vacuum.yh@oerlikon.com
sales.vacuum.yh@oerlikon.com

Oerlikon
Leybold Vacuum
Japan Co., Ltd.
Osaka Sales Office
5-13, Kawagishi-machi
Suita-chi
Osaka 564-0037
Phone: +81-6-6393-5211
Fax: +81-6-6393-5215
info.vacuum.os@oerlikon.com
sales.vacuum.os@oerlikon.com

Oerlikon
Leybold Vacuum
Japan Co., Ltd.
Tsukuba Technical Service Center
Kogyo Danchi
21, Kasuminosato,
Ami-machi, Inashiki-gun
Ibaraki-ken, 300-0315
Service:
Phone: +81-298 89 2841
Fax: +81-298 89 2838
info.vacuum.iik@oerlikon.com
sales.vacuum.iik@oerlikon.com

Korea
Oerlikon
Leybold Vacuum Korea Ltd.
3F, Jellzone 2 Tower, 159-4
Jeongja-Dong, Bundang-Gu
Sungnam-Si, Gyeonggi-Do
Korea 463-384
Sales:
Phone: +82-31 785 1367
Fax: +82-31 785 1359

623-7, Upsung-Dong
Cheonan-Si
Chungcheongnam-Do
Korea 330-290
Service:
Phone: +82-41 589 3035
Fax: +82-41 588 0166

Singapore
Oerlikon
Leybold Vacuum
Singapore Pte Ltd.
No.1, International Business Park
B1-20B, The Synergy
Singapore 609917

Sales and Service:
Phone: +65-6303 7000
Fax: +65-67730 039
sales.vacuum.sg@oerlikon.com
service.vacuum.sg@oerlikon.com

Taiwan
Oerlikon
Leybold Vacuum Taiwan Ltd.
No 416-1, Sec. 3
Chung-Hsin Rd., Chu-Tung
Hsin-Chu, Taiwan, R.O.C.
Sales and Service:
Phone: +886-3-500 1688
Fax: +886-3-583 3999
sales.vacuum.hc@oerlikon.com
service.vacuum.hc@oerlikon.com

oerlikon
leybold vacuum

www.oerlikon.com